

Diamond Microdermabrasion Ultrasound Facial Machine

SE-EM031

User Manual



Accessories Introduction



pores, aging spots, pigment irregularities,
stretch marks, sun damage, wrinkles, etc




reduce fine lines, wrinkles, sagging, clear up
acne, dissolve cellulite, shrink pores, tighten
skin, improve orange peel, etc



improve facial blood circulation, release pain,
smooth wrinkle, tighten pores, improve
orange peel, etc

Control Panel

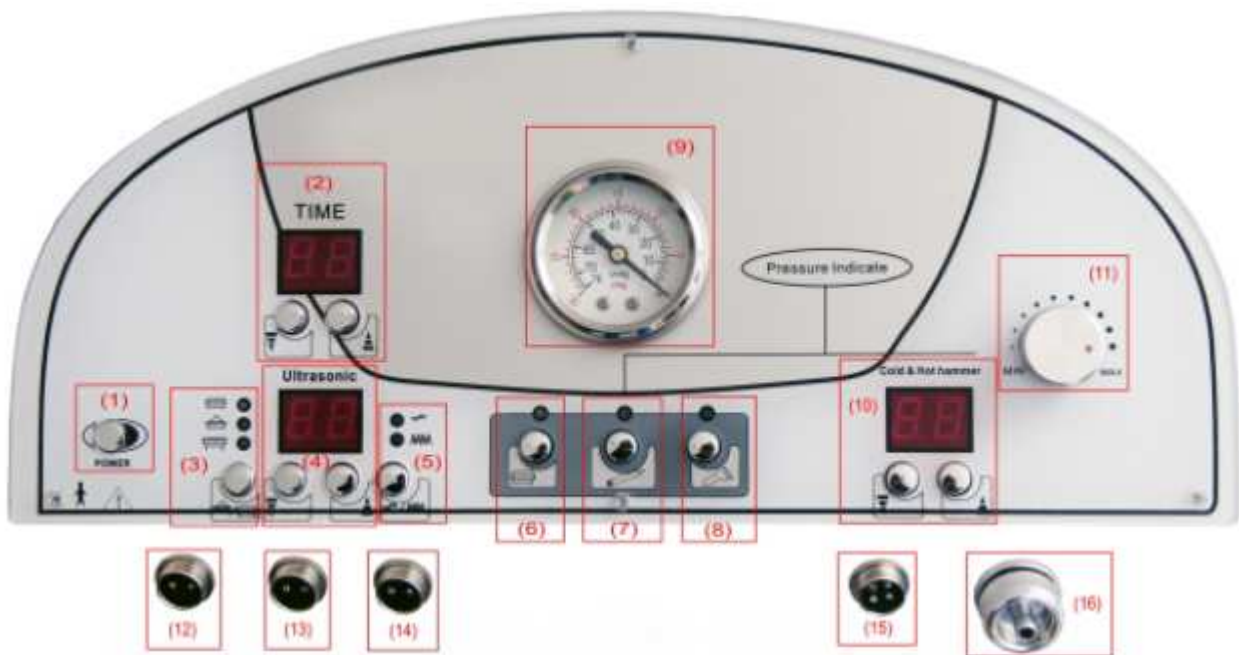


- | | |
|-------------------------------------|--------------------------------------------|
| (1) Power Switch | (9) Vacuum Intensity Display |
| (2) Time Display | (10) Hot & Cold Massager Intensity Display |
| (3) Ultrasound Heads Selection | (11) Vacuum Intensity Adjustment |
| (4) Ultrasound Intensity Adjustment | (12) Body Ultrasound Probe Connection |
| (5) Ultrasonic Waveform Selection | (13) Eye Ultrasound Probe Connection |
| (6) Ultrasound Selection | (14) Face Ultrasound Probe Connection |
| (7) Microdermabrasion Selection | (15) Hot & Cold Massager Connection |
| (8) Hot & Cold Massager Selection | (16) Microdermabrasion Hose Connection |

Operation Control

To use the machine, please firstly connect the power cord to the power supply and then turn on the switch at the back. You will hear a short “beep” which means the machine is in powered. Press the Button (1) on the control panel to turn it on and press the buttons in (2) to select a treatment time. Now the machine is ready for the function you are going to select.

-Diamond Tips Microdermabrasion

Choose a suitable diamond tip and screw on the dermabrasion wand with a piece of cotton filter padding in, connect it to the vacuum hose and connect the other end of the hose to the vacuum inlet at (16) on the panel. Press button (7) to start the vacuum pump and turn the knob at (11) to adjust the vacuum strength. Now it is ready to use on the target area. If you want to pause or stop, press button (7) again.

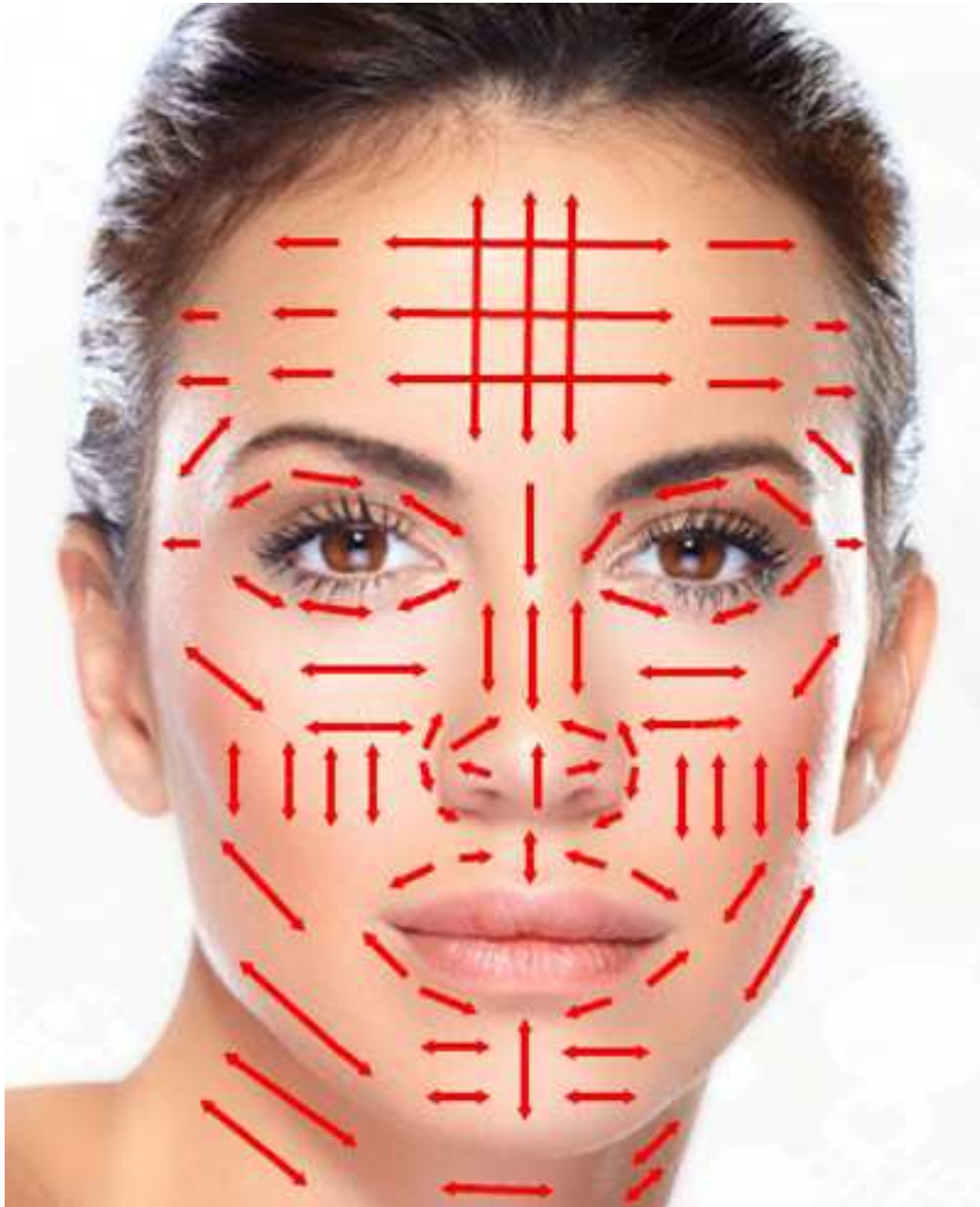

Treatment Schedule Suggestion

- Roughness Peel: 6-10days/treatment, total 3-5 treatments
- Fine Lines: 7-13days/treatment, total 6-8 treatments
- Sagging Neck: 12-14days/treatment, total 3-5 treatments
- Acne Scars: 7-12days/treatment, total 8-10 treatments
- Hyperpigmentation: 7-12days/treatment, total 8-10 treatments

-Ultrasound Facial Treatment

Connect the cords on the ultrasound facial probes and screw the other end of the cord on the ultrasound head inlets on the panel. Inlet (12) is connected to body ultrasound probe, Inlet (13) is connected to eye ultrasound probe, Inlet (14) is connected to face ultrasound probe. Press button (6) to switch to ultrasound function, next press button (3) to select the suitable probe for the target area, then select the intensity in (4) and the waveform in (5). Now the probe is ready to use on the target area. To pause or stop ultrasound treatment, press button (6) again.



-Hot & Cold Massage Hammer

Connect the hot&cold massager to Inlet (15) on the panel and press button (8). Then the fan in the hot&cold massager is running. Turn on the switch on the hot&cold massager so that both the warm side and cold side can start to work. Adjust the intensity at (10) and wait for 1 minute(for cold side) to 2 minute(for warm side) and then it is ready to use on the target area. Turn off the switch on the hot&cold massager if you want to pause. Press button (8) again if you want to stop the treatment.

Technical Parameters

- Voltage: 110V/60Hz, 220V/50Hz
- Power: 135 W
- Packaging Size: 65 x 60 x 32 cm
- G.W.: 10 kg